


SAMPLIFY Application



Integrated liquid level detection

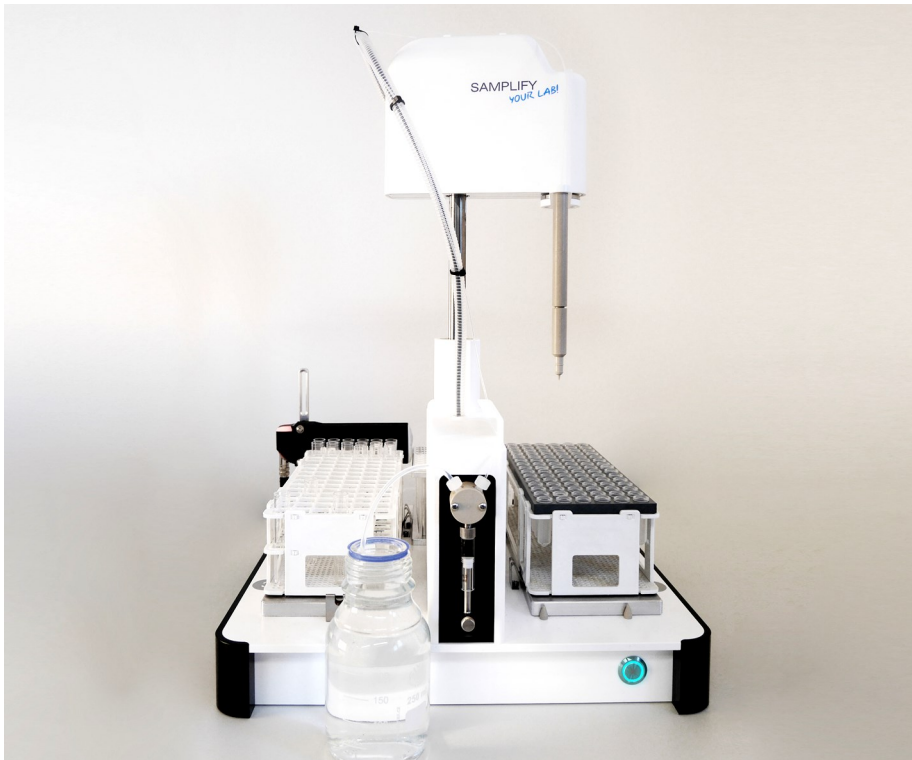

Applications

Sample Preparation for ICP, AAS

No metal Parts in touch with the sample

Working as an automated hand hold pipette with optical control of the required volume

Sample Preparation for medical laboratories

No carryover by Usage of Single-use Tips

Gapless retreat-ability

serum preparation, depending on the volume

Features & Benefits

Pipetting of low, unknown or changing sample volumes, or processes requiring a high degree of control and retreat ability.


Sample preparation depending on the sample volume. Tracking via Barcode with transfer of 60 barcodes by one scan.

Predestined configuration for working with low and varying sample volumes.

Guaranteed no carry over by usage of Single-use-Pipette Tips (Adjustable on all common pipette tips by providing individual adapters)

Level Control for varying sample volumes

- Pipetting the required volume
- Checking the volume
- If the volume is correct: diluting.
If not: pipetting a lower volume
- Checking the lower volume
- If the volume is correct: diluting with a lower liquid amount
If no: re-pipetting the sample and marking the vial in the software

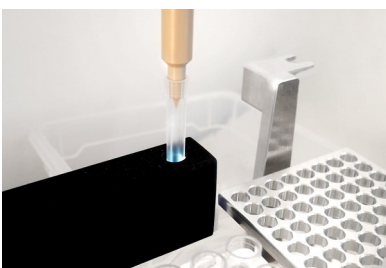

nevoLAB GmbH
Am Gehrenbach 8
D 88167 Maierhöfen
phone: +49 (0) 8383 929 566 0
fax: +49 (0) 8383 929 566 30
email: info@nevolab.de
www.nevolab.de